

# **XRF-XFS-XAS- Auto v.1.0 (beta) Manual**

ALS XFM Beamline 10.3.2

Sirine C. Fakra (April 2024)



## XRF-XFS-XAS- Auto



XRF-XFS-XAS-Auto Copyright (c) 2024

This automated data analysis software consists of three parts for

1.  $\mu$ X-ray fluorescence mapping (in blue)
2.  $\mu$ X-ray fluorescence spectroscopy (in green)
3.  $\mu$ X-ray absorption spectroscopy (in pink)

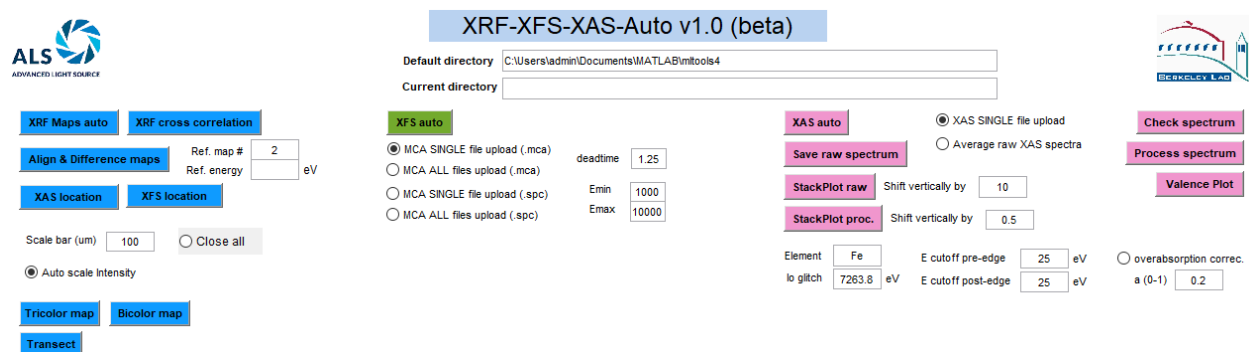


Fig.1

The analysis features for each part are described below.

## 1. Micro- X-ray Fluorescence mapping ( $\mu$ XRF) processing

Only maps that have been deadtime-corrected can be processed.

### 1.1 Elemental XRF maps

Maps are read and then normalized by the dwell time and  $I_0$  (incident beam intensity). To display the elemental maps start by clicking on the button “XRF maps auto” (**Fig. 2**).


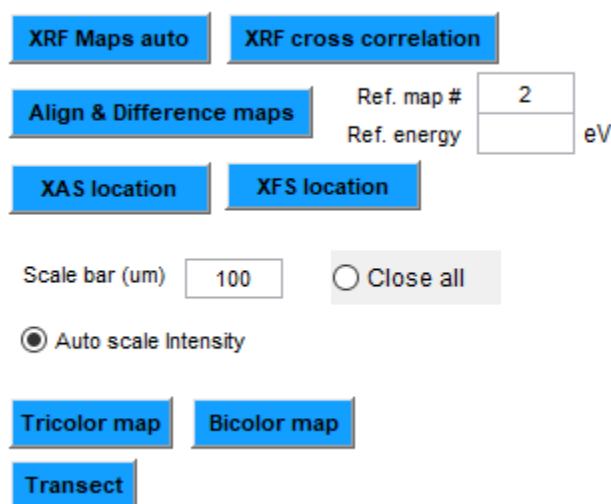



Fig. 2

Select the .xrf map file that you want to open (**Fig. 3**). You can open up to 2 map files maximum.

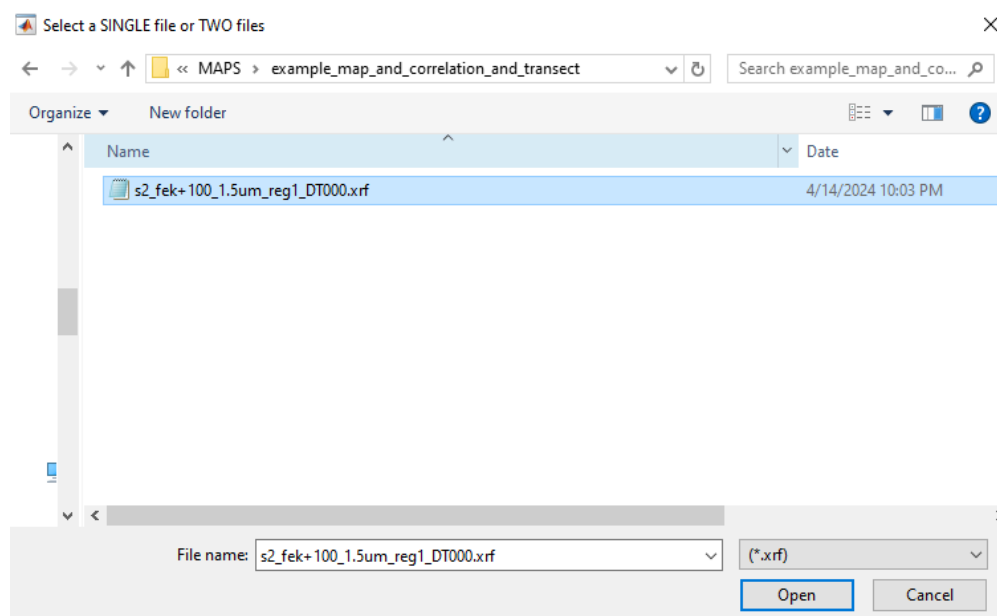


**Fig. 3**

This window will then pop asking which elemental map you want to display (**Fig. 4**), in this case Fe, Ca and K.

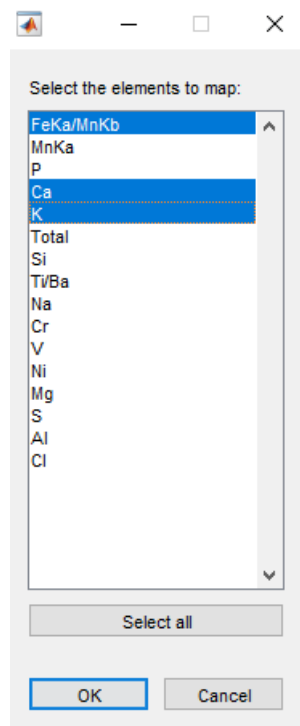

**Fig.4**

Elemental maps are then displayed (**Fig.5**) and automatically saved as .bmp files within the same folder.

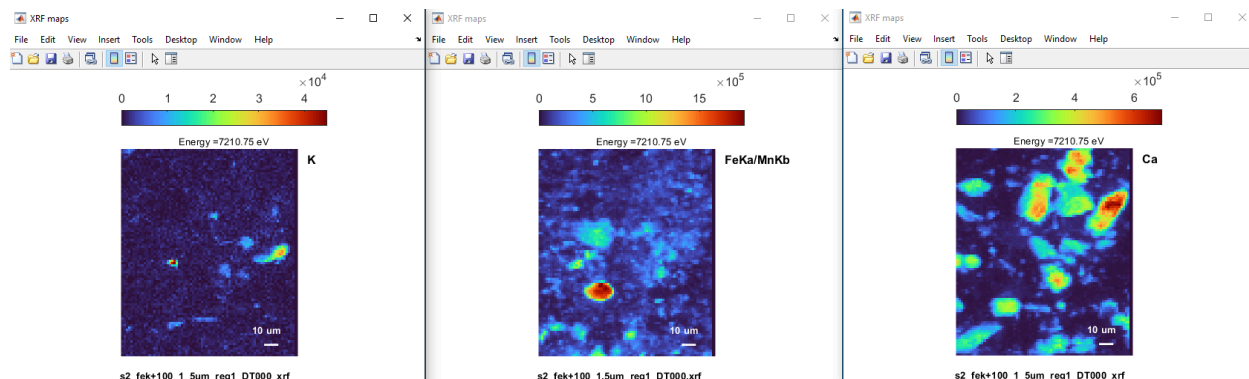

**Fig. 5**


In case you need to change the scale bar, just type another value next to “scale bar” (**Fig.2**) and click on “XRF map auto” again, the current maps will close automatically, the new ones will be displayed and the previous .bmp maps will be overwritten. Once you’re done, click on the “close all” button (**Fig.2**) to close all maps.

## 1.2 Tricolor and bicolor maps

Tricolor and bicolor maps use the RGB convention: red=1st element, green=2nd element and Blue=3rd element.

### 1.2.1 Tricolor map


First click on the “XRF map auto” button (**Fig.2**) and follow the steps in section 1.1. You must have only three elemental maps opened, in this case Fe, Ca and K, so close all maps that you don’t need. Press the “tricolor map” button, a composite map of these three elemental maps will be displayed (**Fig. 6**), with iron in red, calcium in green and potassium in blue. The tricolor map is saved automatically within the same folder.

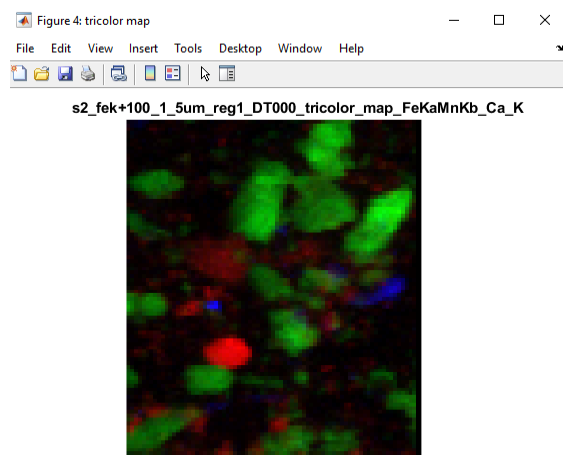

**Fig. 6**

### 1.2.2 Bicolor map

First click on “XRF map auto” (**Fig.2**) and follow steps in section 1.1. You must have only two elemental maps opened, in this case Fe and Ca, so close all maps not needed. Press the “bicolor map” button, a composite map of these two elemental maps will be displayed (**Fig.7**), with iron in red and calcium in green. The bicolor map is saved automatically within the same folder.

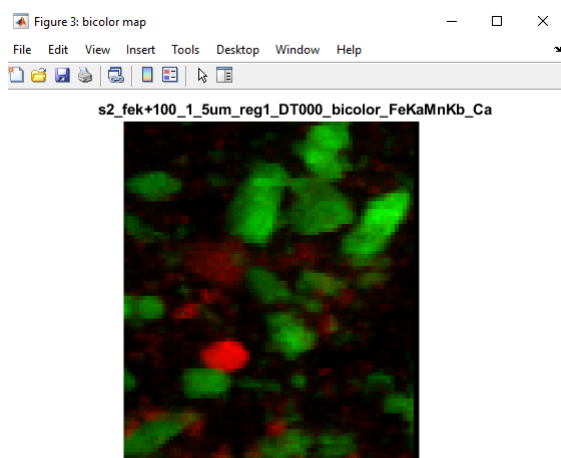


**Fig. 7**

### 1.3 XRF cross-correlation


Click on “XRF map auto” (**Fig.2**) first, pick a map file, you can open up to 2 map files. To check whether two elements are correlated or not, in this case Fe and Ca (**Fig. 7**), press on the “XRF cross-correlation” button.

A message window will pop-up (**Fig. 8**). If you want to check the correlation between 2 elements from the same map, click “No” as illustrated in this example, otherwise press “Yes” if you want to check the correlation between 2 elements coming from 2 different maps, typically recorded at different energies.

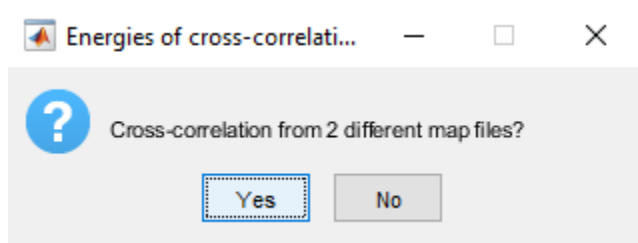


**Fig. 8**

Then select the first element (**Fig. 9**), press “OK” and select the 2nd element, and press “OK”.

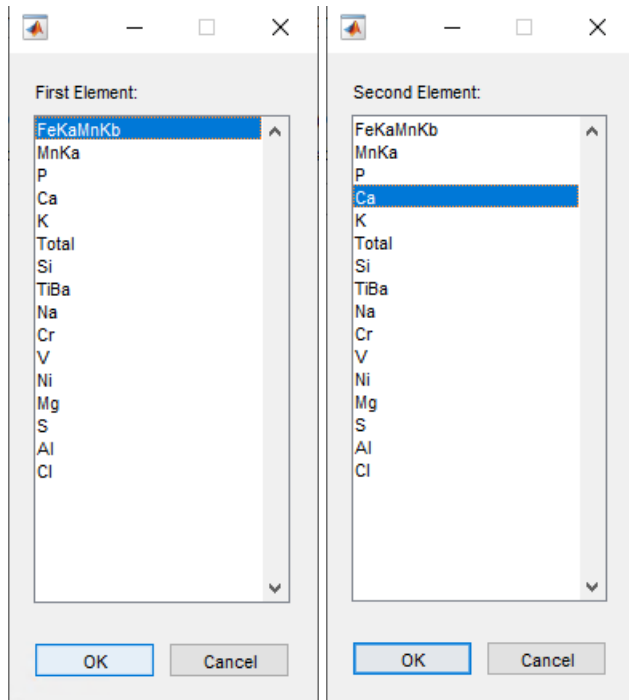

**Fig. 9**

The correlation plot between these 2 elements is displayed (**Fig. 10**). The plot is saved automatically as CorrelationPlot\_element1\_element2 within the same folder. To check correlation between other elements, press on the “XRF cross-correlation” button again and repeat the steps above. You can generate as many correlation plots as you want. Press “Close all” when you’re done.

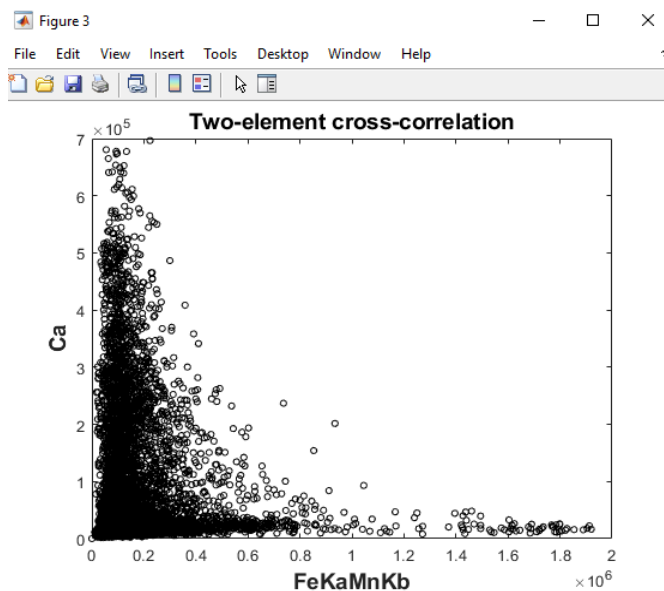

**Fig. 10**


## 1.4 Align and Difference maps

To obtain “difference” maps, you must first align all maps to a reference map in the set, typically one recorded below the edge. For this, in the reference map number box (**Fig. 11**), type the number of the reference map, with scan 000 corresponding to map# 1. Click on the “Align and Difference maps” button.




**Fig. 11**

Select any map in the folder that contains all the maps you want to align (**Fig. 12**) then click “open”

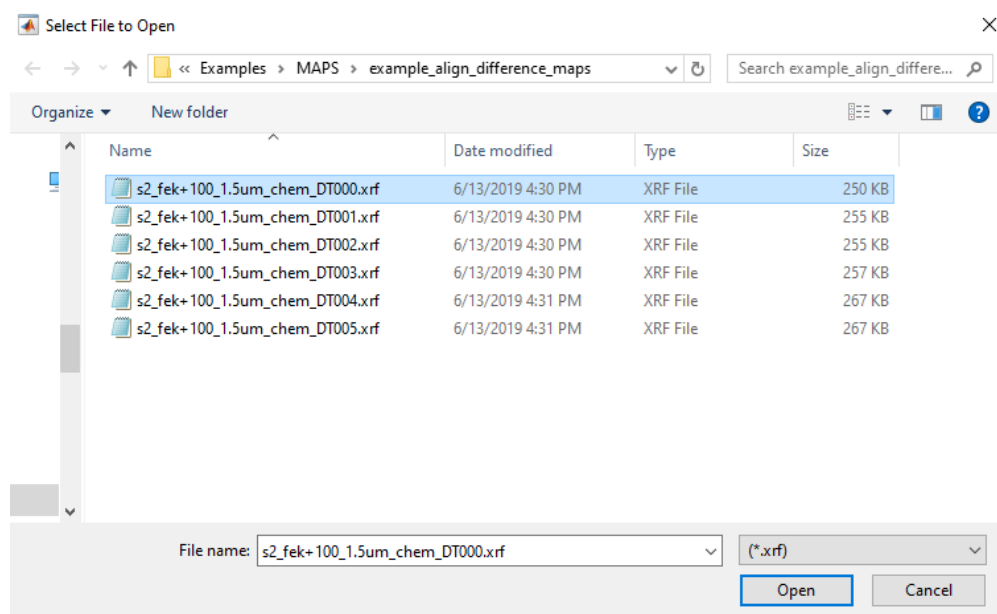


**Fig. 12**

In the “Pick elemental map for the offset” window (**Fig. 13**), choose an elemental map, K in this example, where you have small features, for instance a particle, that you can use to align the maps. Once the K map is displayed, zoom on a small and bright feature, then press “ok” to continue. The location of the feature is used to properly register the maps at different energies. If you encounter the message “too low contrast for accurate alignment”, try another feature on the same map where the contrast is higher or try another element for the offset map. Then, pick the element for the “difference” map, Fe in this example, and click “OK”.

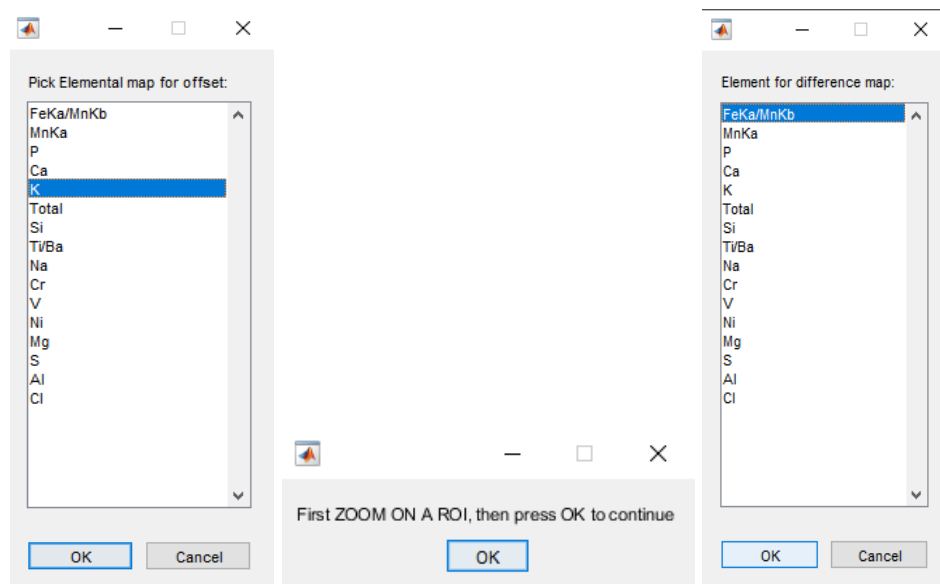

**Fig. 13**

Difference maps are then displayed (**Fig. 14**) and automatically saved within the same folder.

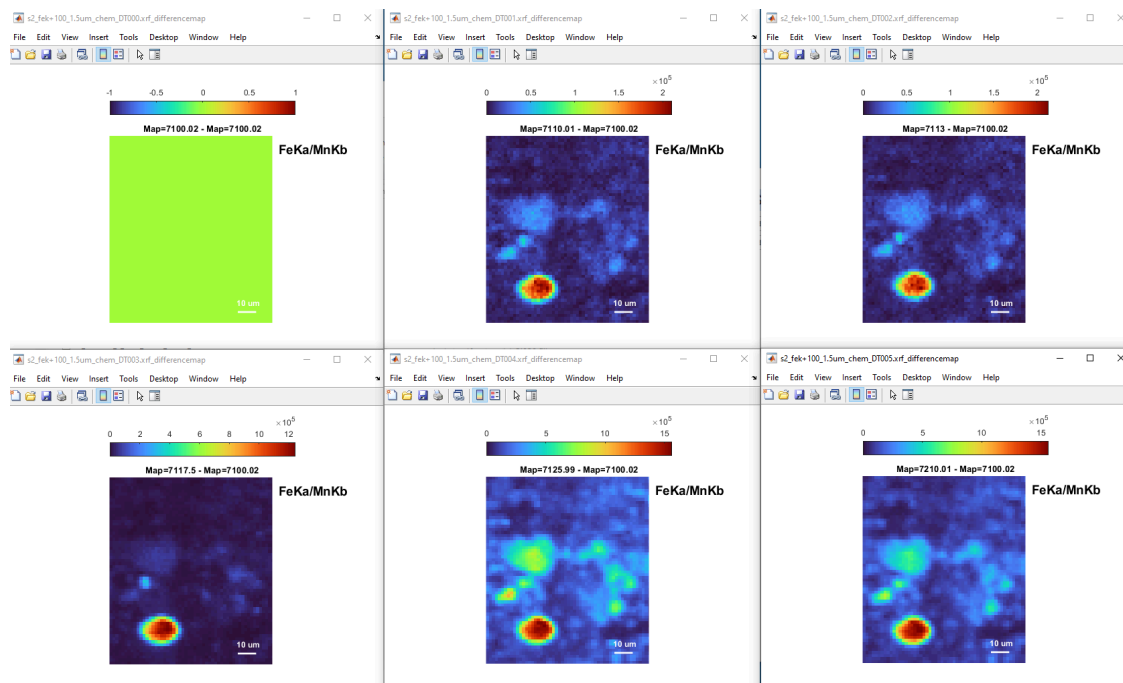


**Fig. 14**

### 1.5 Transect


To check the profile along a transect in the region of your sample, click on “XRF maps auto” to display elemental maps, close the maps that you don’t need, then click on the “transect” button (**Fig.2**). A message window will pop-up (**Fig. 15**), press “OK” to activate the cursor on the map



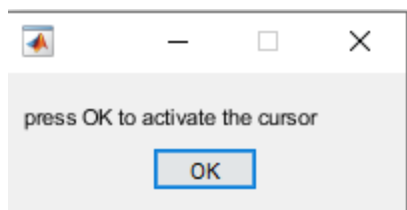

**Fig. 15**

Then click on a starting point and an ending point of the transect on the map. Be sure to choose points that are inside the map area only (**Fig. 16**).

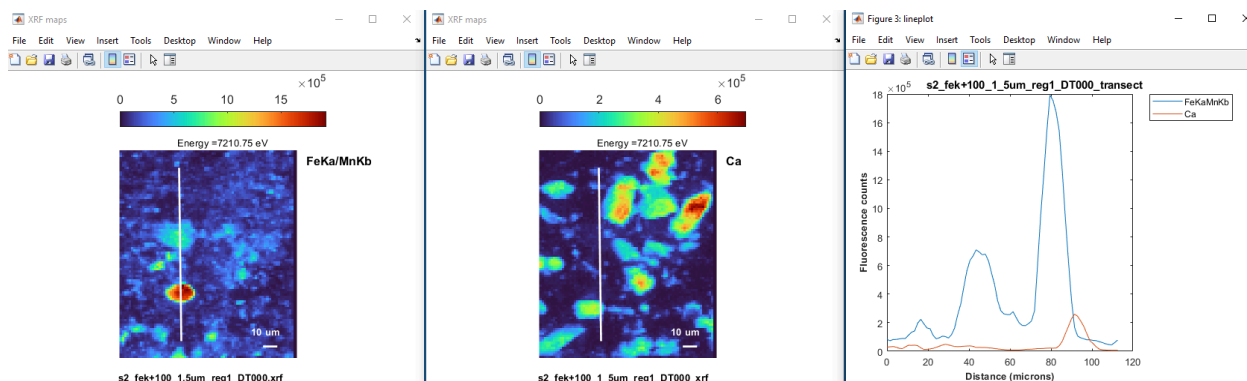

**Fig. 16**


The transect plot and annotated maps with transect location are automatically saved as “filename\_transect#.bmp” in a folder named “filename\_transect#”.

To do another transect, just press the “transect” button again, and repeat steps above.

Transects are obtained from all elemental maps displayed. If there’s an elemental transect that you don’t want, just close the corresponding map and repeat. Once you’re done, press the “close all” button.

## 1.6 QXAS location


To display the locations where you collected XAS spectra, click on “XRF map auto” to display your maps, then press on the “XAS location” button (**Fig. 2**). Select a QXAS file (.qx), the QXAS location will be shown on the maps as a white-empty circle, labeled with the scan# above it (**Fig. 17**). These maps are automatically saved as .bmp in a folder named “filename\_XAS”.

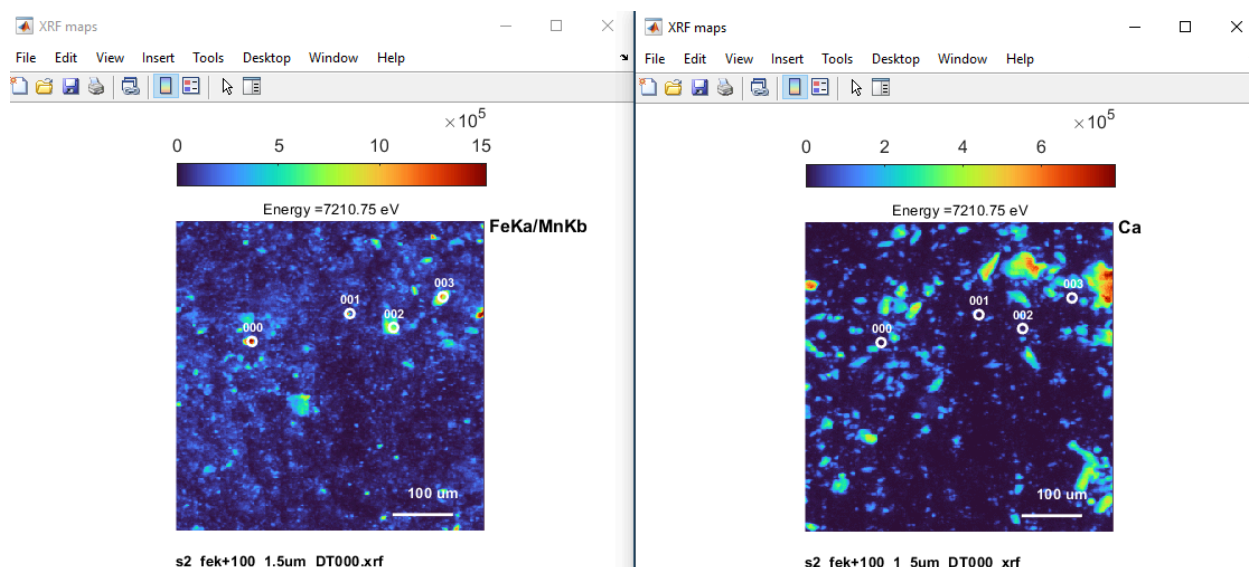

**Fig. 17**

### 1.7 XFS location

To display the locations where you collected XFS spectra, click on “XRF map auto” to display your maps, then press on the “XFS location” button (**Fig.2**). Select a MCA file (.mca), the XFS location will be shown on the maps as a white-empty circle, labeled with the scan# above it (**Fig. 18**). These maps are automatically saved as .bmp in a folder named “filename\_XFS”.

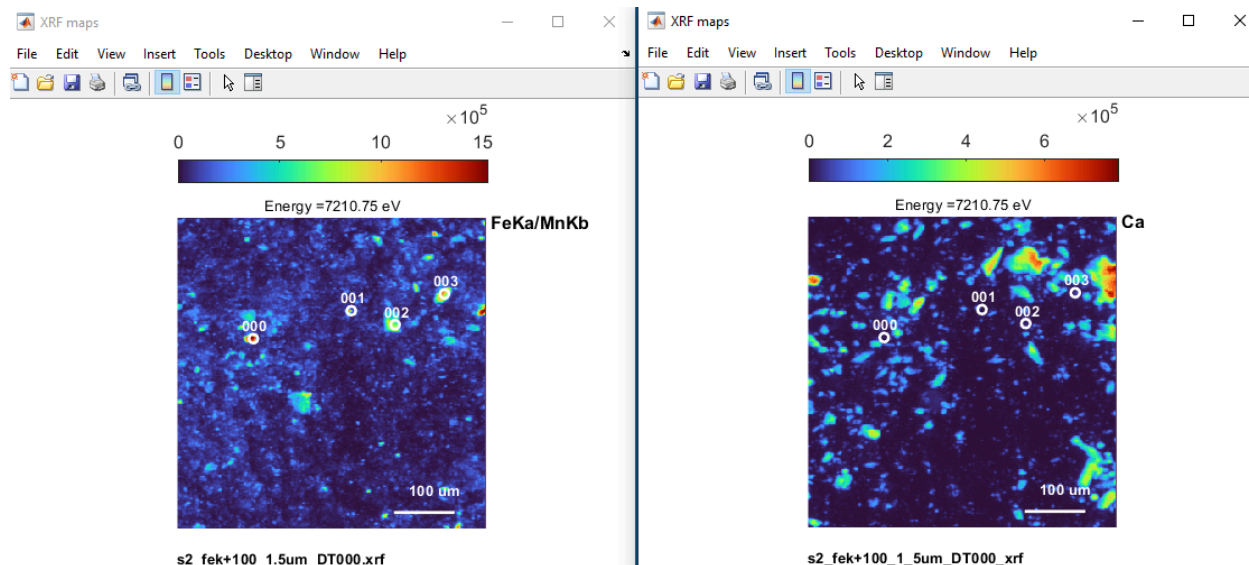

**Fig. 18**

## 2. Micro- X-ray fluorescence ( $\mu$ XFS) processing

You have the option of either processing raw (.mca) spectra, or displaying processed spectra (.spc). The first option is preferred, as automatic energy calibration is performed in this case.

**XFS auto**

☒ MCA SINGLE file upload (.mca)  
☐ MCA ALL files upload (.mca)  
☐ MCA SINGLE file upload (.spc)  
☐ MCA ALL files upload (.spc)

deadtime   
Emin   
Emax

Fig. 19

Fluorescence emission lines corresponding to K-L3 and L3-M5 transitions are labeled on all plots. Ka emissions lines are labeled, but not the Kb to avoid too many labels displayed. The elastic/compton peak is not labeled.

### 2.1 Single raw file (.mca)

To process raw files (.mca), enter the deadtime correction value, in the “deadtime” window (Fig. 19), then choose the “MCA single file upload (.mca)”. The energy is calibrated using two emission lines ArKa and FeKa in this example, with  $E_{\text{calib}}(\text{Ar}) = 2957.7$  eV,  $E_{\text{calib}}(\text{Fe}) = 6403.84$  eV. Fig. 20 illustrates the energy calibration.

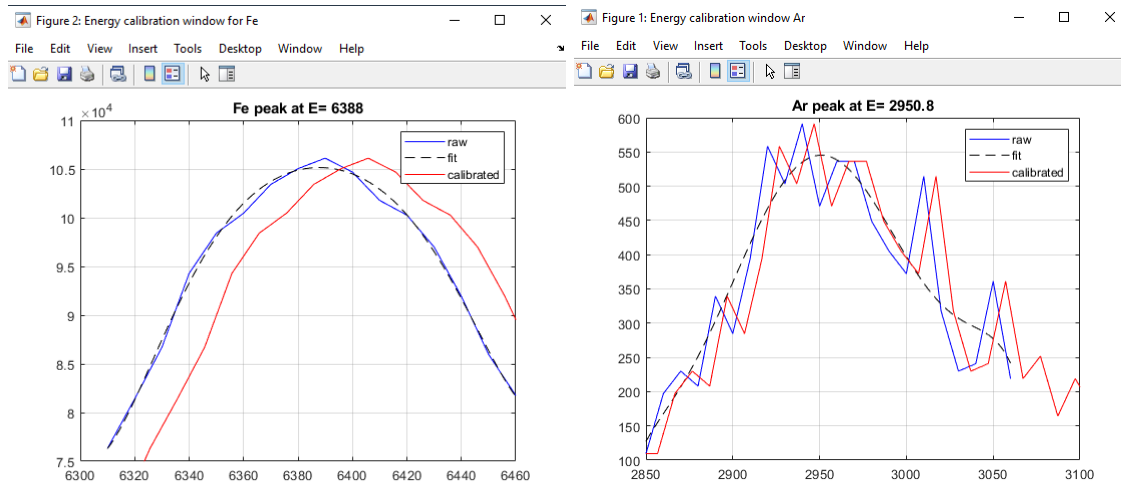


Fig. 20

The spectrum is then displayed in base 10 logarithmic scale (Fig. 21), with peaks labeled (Ka and La lines only).

This plot is saved as .bmp within the same folder.

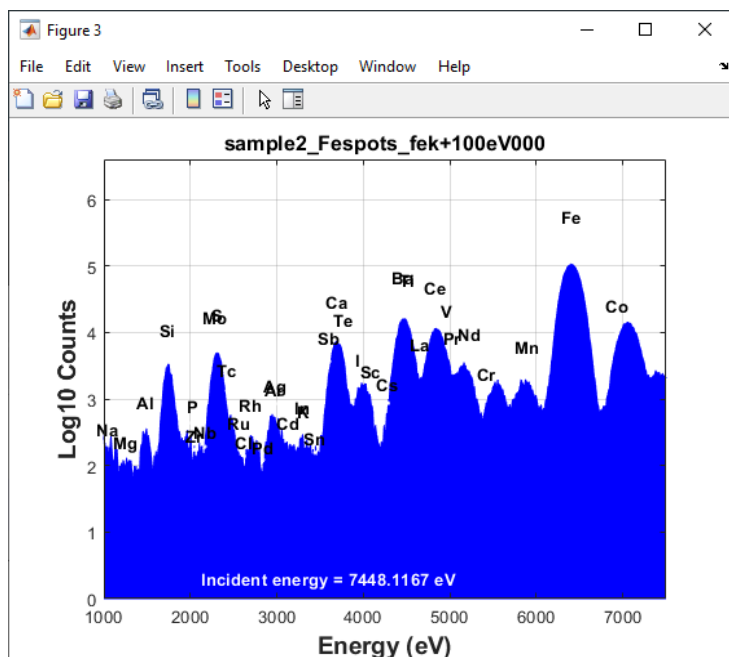


Fig. 21

## 2.2 Single processed file (.spc)

The .spc file comes deatime-corrected, acquisition time and Io normalized. Choose the “MCA single file upload (.spc)” option (**Fig. 19**), select a .spc file and the XFS spectrum will be displayed in log10 scale with fluorescence emission lines labeled (K $\alpha$  and L $\alpha$  lines), see **Fig. 22**. This plot is automatically saved as .bmp within the same folder.

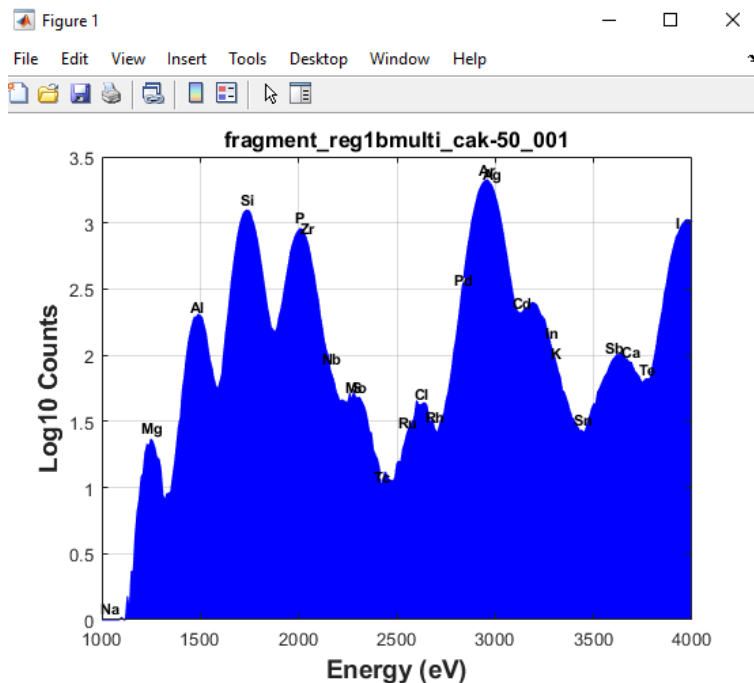


Fig. 22

### 2.3 Multiple raw (.mca) and processed (.spc) files

If you want to analyze multiple raw spectra (.mca files), click on the “MCA all files upload” button (**Fig.19**) then choose any .mca file in the folder where all your files are. Same for the .spc files, except choose the “SPC all files upload”.

All files contained in the folder will be plotted and the plots saved automatically as .bmp within the same folder.

### 3. Micro- X-ray absorption spectroscopy ( $\mu$ XAS)

This software is intended for the analysis of  $\mu$ XANES and extended- $\mu$ XANES spectra only. You can open both .qx files and .dat files. Note that incomplete QXAS spectra cannot be processed.

The screenshot shows a software interface for X-ray absorption spectroscopy analysis. It features several buttons and input fields. At the top left, there is a button labeled "XAS auto". Below it are buttons for "Save raw spectrum", "StackPlot raw", and "StackPlot proc.". To the right of these buttons are two radio buttons: "XAS SINGLE file upload" (selected) and "Average raw XAS spectra". Further right are buttons for "Check spectrum", "Process spectrum", and "Valence Plot". Below the buttons, there are input fields for "Shift vertically by" with values "10" and "0.5". At the bottom, there are input fields for "Element" (Fe), "E cutoff pre-edge" (25 eV), "E cutoff post-edge" (25 eV), "Io glitch" (7263.8 eV), and "overabsorption correc. a (0-1)" (0.2).


XAS auto	<input checked="" type="radio"/> XAS SINGLE file upload	Check spectrum
Save raw spectrum	<input type="radio"/> Average raw XAS spectra	Process spectrum
StackPlot raw	Shift vertically by 10	Valence Plot
StackPlot proc.	Shift vertically by 0.5	
Element Fe	E cutoff pre-edge 25 eV	<input type="radio"/> overabsorption correc.
Io glitch 7263.8 eV	E cutoff post-edge 25 eV	a (0-1) 0.2

**Fig. 23**

### 3.1 Displaying and checking raw spectra

#### 3.1.1 Displaying raw spectra

To open a raw spectrum, first click on “XAS SINGLE file upload”, put the element, then press on “XAS auto” (**Fig. 23**). The raw spectrum will be displayed on the “Raw” panel. In this example (**Fig.24**), the absorption spectrum of an iron foil, recorded in transmission mode (one active scaler, #14). Scaler #0 corresponds to the incident beam intensity,  $I_0$ .

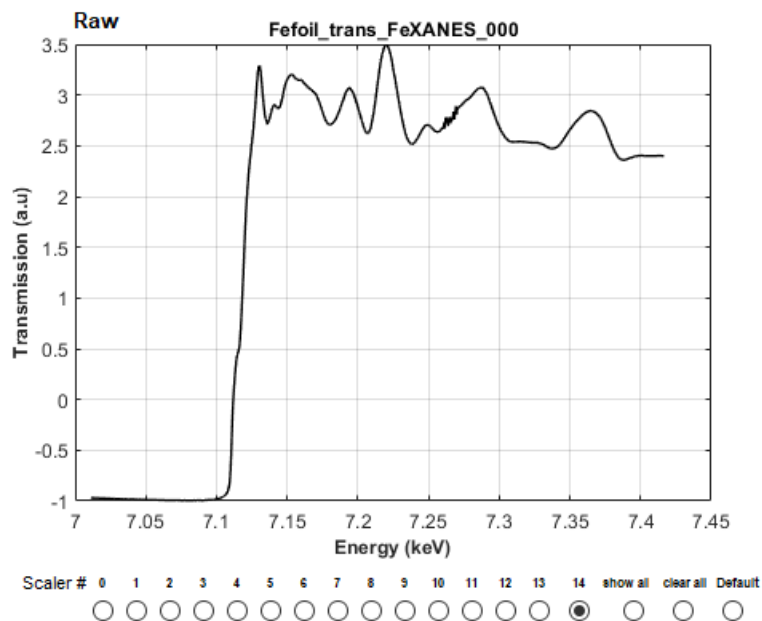

**Fig. 24**

You can select all scalers, unselect them or clear them all by clicking the corresponding buttons. If your data were collected in fluorescence mode with the multi-element SDD detector, several scalers will be active and you can check each channel individually to see if you have any outliers. You can save the raw spectrum by pressing “Save raw spectrum”.

### 3.1.2 Average of raw spectra


To average several raw spectra, recorded consecutively at the same location on an XRF map, select “Average raw spectra”, then press “XAS auto” (**Fig. 23**). Select all spectra you want to average then click “open”. The average will be displayed on the “Raw” panel.

You can save the average by pressing “save raw spectrum”, the .bmp file will be saved in a separate folder “filename\_scan#\_XAS\_average”.

### 3.1.3 Stack Plot raw spectra

This feature is used to quickly check which raw spectrum you want to process further. Press the “StackPlot raw” button, and select all spectra that you want to display (**Fig. 25**). You can change the vertical shift value.


This plot is automatically saved within the same folder as Stack\_Plot\_raw\_#.

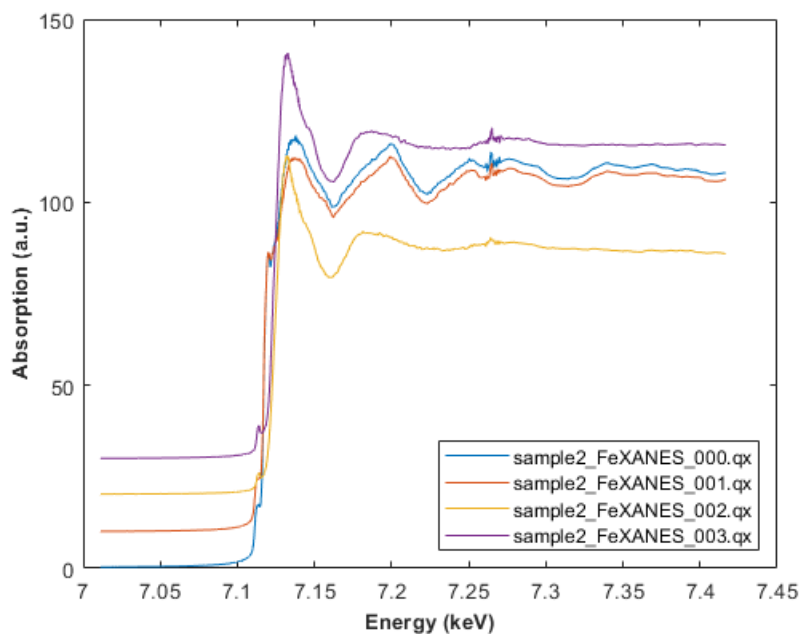

**Fig. 25**

### 3.1.4 Checking QXAS raw spectra

A useful option for checking a QXAS raw spectrum is to look at the shape of the spectrum at the beginning and end of its acquisition. Press “XAS auto” to upload a raw spectrum, then press “Check spectrum”. The average over the first 3 lines and over the last 3 lines of the spectrum containing 20 lines total, will be overlaid, as illustrated in **Fig. 26**. This feature is useful to check for radiation beam-induced damage or for sample drifting for instance.

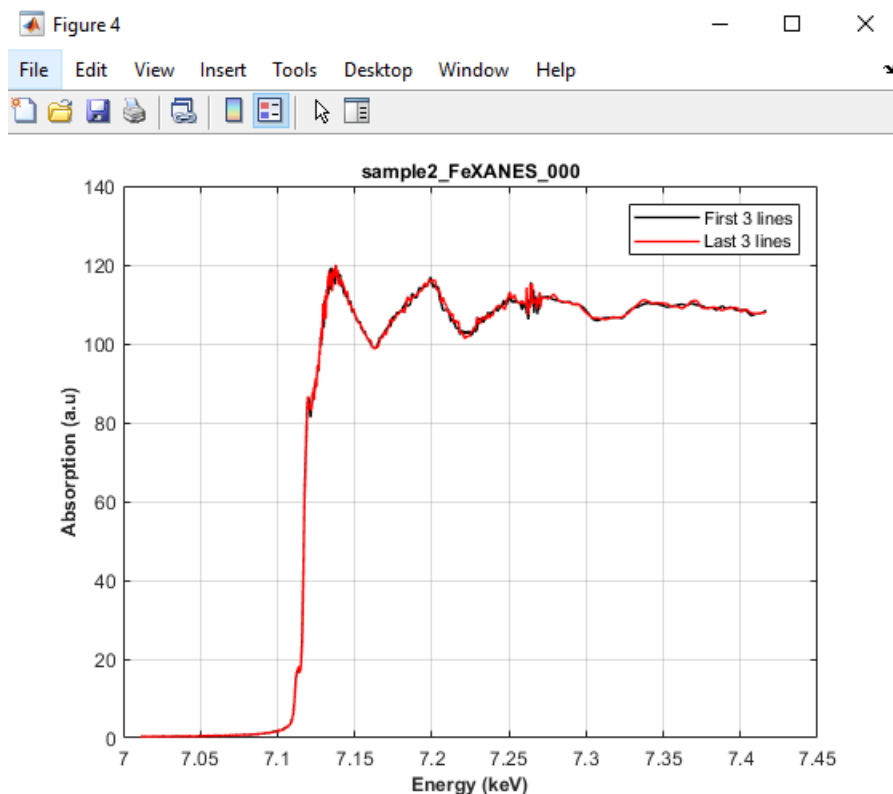

**Fig. 26**

## 3.2 XAS spectrum processing


To process a raw spectrum, make sure to click back on “Default” on the scalers in the Raw panel.

### 3.2.1 Energy calibration

Type the value of the Io glitch in the “Io glitch” window and put the element (**Fig. 23**). Then press “process spectrum”, the spectrum will be automatically calibrated.

### 3.2.2 Pre-edge background subtraction, post-edge normalization

Click on the “process spectrum” button, then three options for the pre-edge fit will pop (**Fig. 27**).

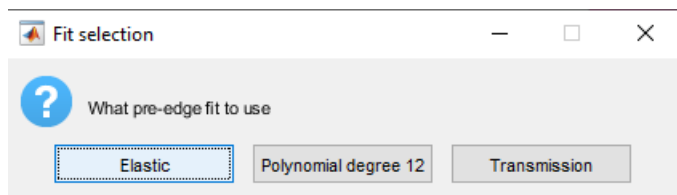

**Fig. 27**

“Elastic”: the default option for spectra recorded in fluorescence mode. Following the equation:



$Y = a + b/(E_0 - E) + c/(E - E_f)^2$ ,  $E_0$ = absorption energy,  $E_f$ = fluorescence emission energy.

Polynomial degree 12: for spectra recorded in fluorescence mode

Following the equation:  $Y = a_0 + a_1 E + a_2 E^2 + a_3 E^3 + \dots + a_{12} E^{12}$

Transmission: for spectra recorded in transmission mode.

Following the equation:  $Y = a E^{-2.7} + b$ ;  $E$ =energy (eV)


In the case of the “Elastic” and “Transmission” options, the program will automatically compute

the  $E$  pre-edge cutoff that minimizes  $NSS2$ , with  $NSS2 = \sum (\mu_{exp} - \mu_{fit})^2 / \sum (\mu_{exp})^2$

Zoom in to inspect the pre-edge, and if not satisfied, you can play with the  $E$  cut off pre-edge value as illustrated in **Fig.28**,  $E$  cut off pre-edge= 35 and  $E$  cut off post-edge= 15.


For the post-edge, a polynomial order 2 is used, following the equation:

$Y = a_0 + a_1 E + a_2 E^2$  Make sure to check that this fit goes well in the middle of the “EXAFS” oscillations.

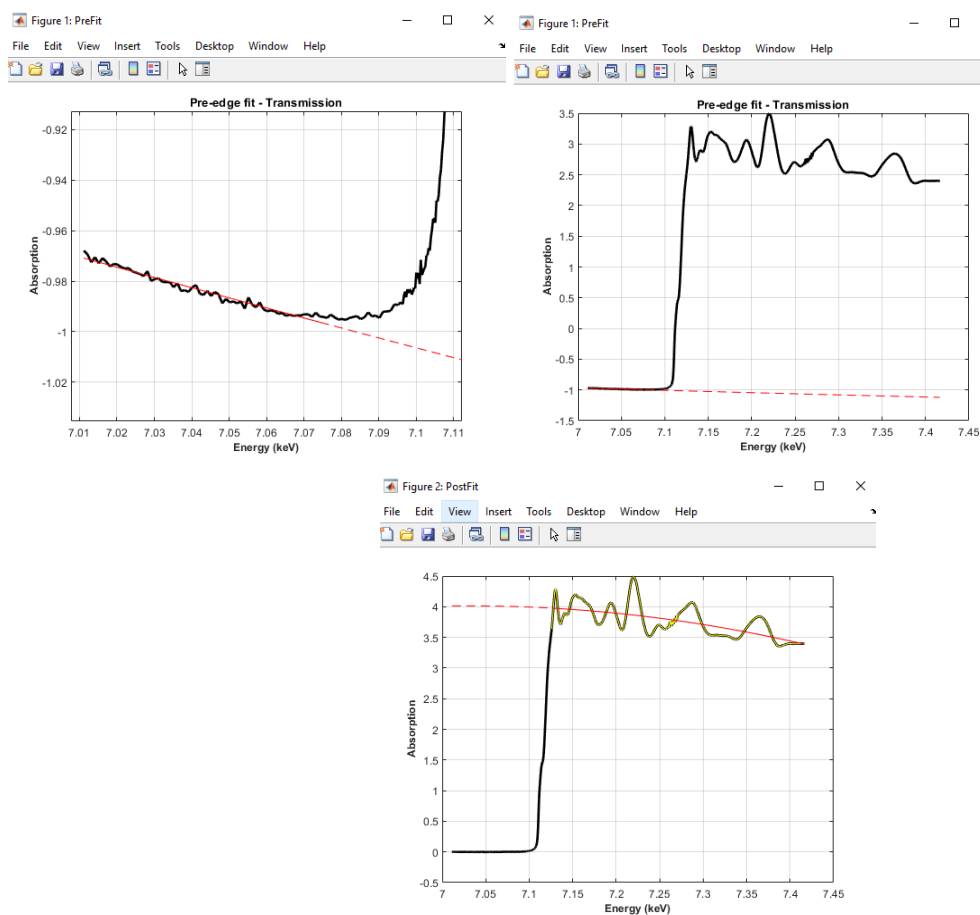


**Fig.28**

The processed spectrum, pre-edge background subtracted and post-edge normalized will be displayed on the “Processed” panel (**Fig. 29**).

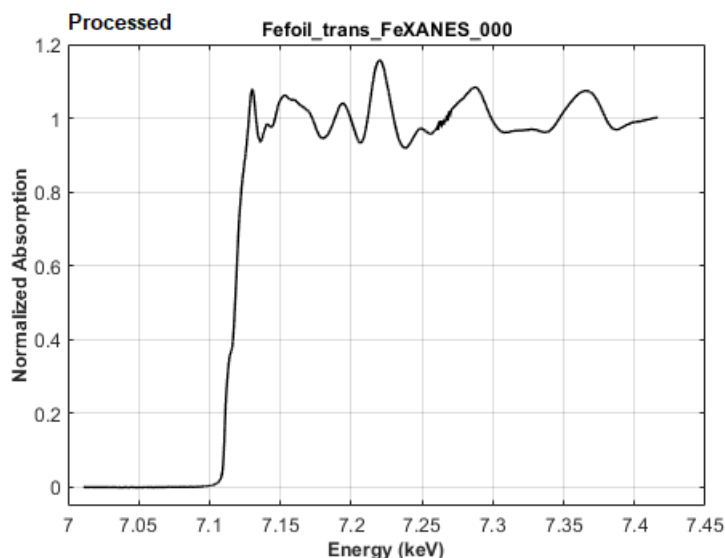


**Fig. 29**

At this stage, a “Select standard folder” dialog window will appear. If you click “cancel”, no fit will be performed and you will then have the option of saving the processed spectrum. Select “Yes” on the following window (**Fig. 30**) to save the spectrum otherwise press “No”.

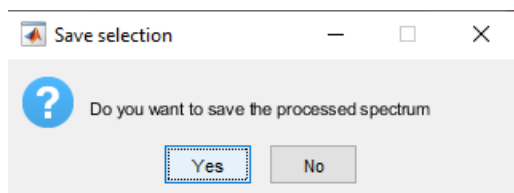


**Fig. 30**

If a Standard folder is selected the program will then automatically proceed with least-square linear combination (LSQ) fitting of the processed experimental spectrum using 1 or more standards and you will be asked to select the folder containing the set of standards as illustrated in **Fig. 31**.

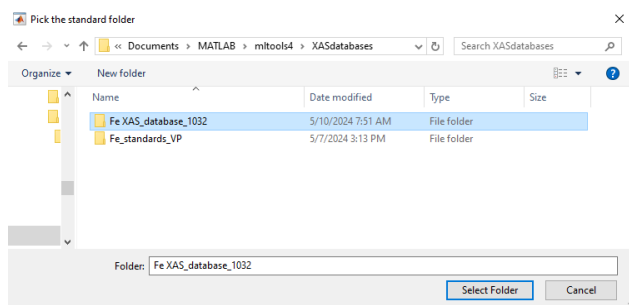


**Fig. 31**

Then, you'll be asked how many standards you want to use for the LSQ fitting. You can use a combination of up to 4 standards (**Fig. 32**).

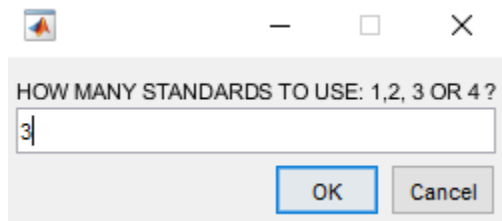

**Fig. 32**

If you're not sure how many standards to use, start with 1, then 2, then 3 and then 4. Look at the sum-sq values, if they differ by less than 10%, then the extra component added is very likely not required. After selecting the number of components to fit, press "ok". The top 5 best LSQ combination fits will be displayed as illustrated in **Fig. 33**, with an iron foil recorded in transmission mode using 1 standard..

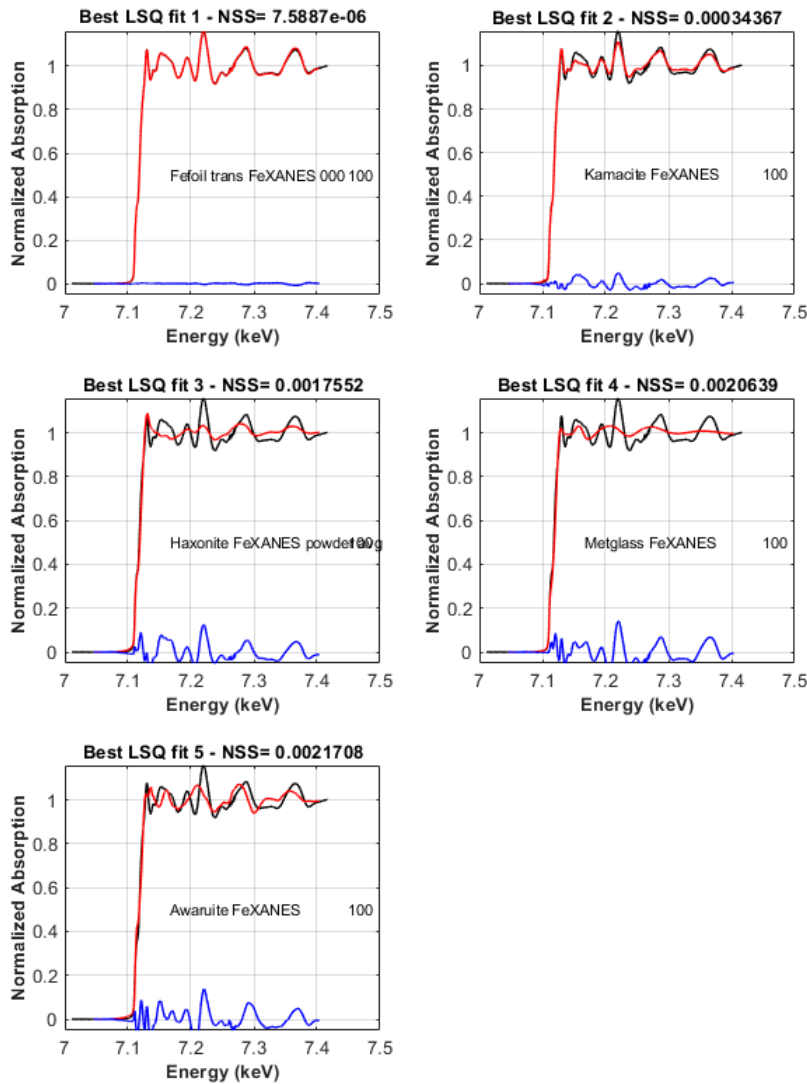

**Fig. 33**


The best mathematical fit is #1 where the NSS1 and NSS2 values are minimum, with

$$NSS2 = \sum (\mu_{exp} - \mu_{fit})^2 / \sum (\mu_{exp})^2$$

$$NSS1 = \sum (\mu_{exp} - \mu_{fit})^2$$

where  $\mu$  = normalized absorption.

Note that when using up to 3 components it takes less than 10 sec for the fits to be displayed while when using 4, it takes about 20 minutes! Once the fits are displayed, always inspect the spectral fine features and look at the NSS values and how they differ and choose your fit. A dialog window will ask you to choose between the top 5 best fits (**Fig. 34**). In this example, we chose Fit#1 which is not only the best mathematical combination but also what makes sense.

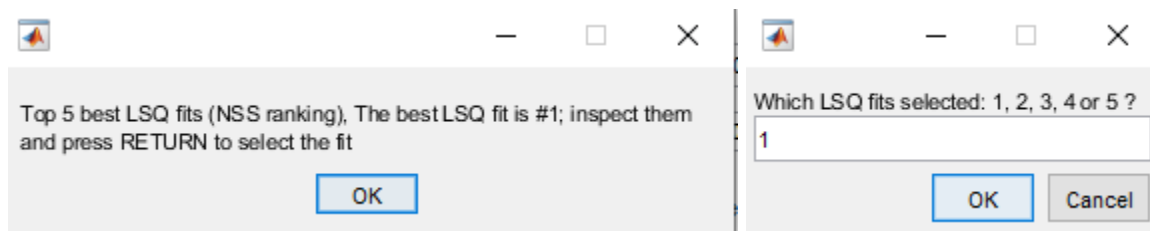


**Fig. 34**

A table showing the major components, their percentage and NSS values will be displayed and saved automatically as an excel and a pdf file (**Fig. 35**) as well as the fits and residuals (**Fig. 36**)

Figure 3: Major Component

Component Name	% in fit	SSN1	SSN2
Fefoil trans FeXANES 000	100	0.0091255	7.5887e-06

**Fig. 35**

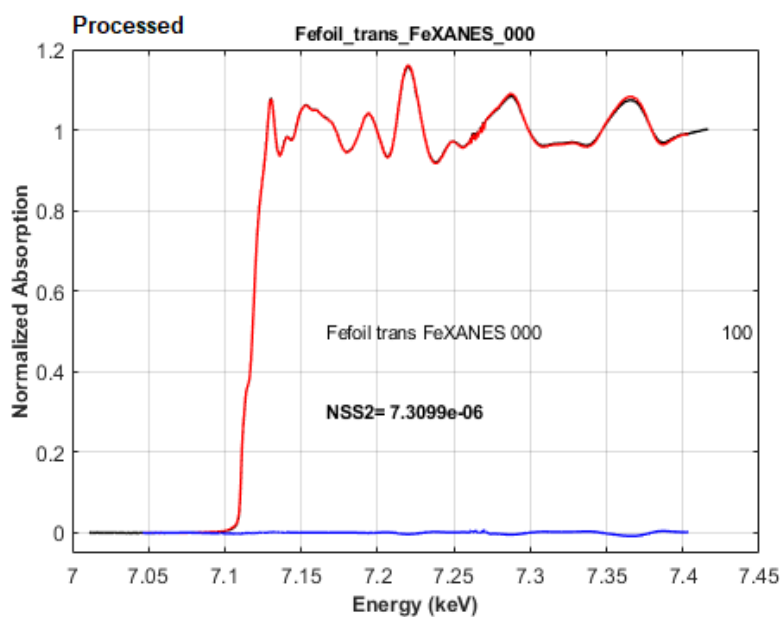

**Fig. 36**

The processed experimental spectrum (in black), the fit (in red) and residual (in blue) will be displayed on the “Processed” panel. Several files are automatically saved in the same folder named “filename\_XAS”:

- Plot in a bmp format
- excel files of the best fit and the processed spectrum

- txt files of the processed spectrum

### 3.5.3 Over-absorption correction

If the fit of the experimental spectrum recorded on a “hot” point is not satisfactory, apply an over-absorption correction factor, “a”, to the spectrum to correct for distortions (**Fig. 37**)



**Fig. 37**

This factor can vary between 0 and 1, however is typically <0.5. The over-absorption correction factor **a** is applied following the equation:  $Y = Y/(a+1-Y*a)$ .

After the over-absorption value is entered, press “Process spectrum” and follow the steps described in the previous section.

### 3.5.3 Valence scatter Plot

The scatter plot is used to quickly classify XANES spectra according to metal valence state. It computes the normalized absorption values,  $\mu(E)$ , at two specific energies  $E_1=7117.5$  eV and  $E_2=7113$  eV. A scatter plot ( $\mu_1, \mu_2$ ) is then obtained, with  $\mu_1=\mu(E_1)$  and  $\mu_2=\mu(E_2)$ , where one datapoint corresponds to one spectrum. This method has been previously described elsewhere (2 references), it has been revised since and expanded to other metals and metalloids.

Press “Valence Plot” to access this function; the window below will appear (**Fig. 38**).

#### XFM10.3.2- XANES-derived Valence Plot

Default directory:  
C:\Users\admin\Documents\MATLAB\mltools4XU\

Select Element: **Fe** (current element)

Open Data Files: Name:  Number of files:

Plot: Undo

Parameters:  
 Mu3:  7110  
 Mu2:  7113  
 Mu1:  7117.5

Plot Parameters:  
☐ Hexagon ☐ Solid ☐ Red ☐ Green ☐ Orange  
☒ Square ☒ Empty ☐ Blue ☐ Cyan  
☐ Circle ☐ Black ☐ Magenta

Font Size:  4


Xmin:  0.025 Xmax:  1 Ymin:  0.01 Ymax:  1

☐ auto axis

CLOSE


**Fig. 38**

Press “select element” and pick one of the three available, for instance here, Fe (**Fig. 39**).

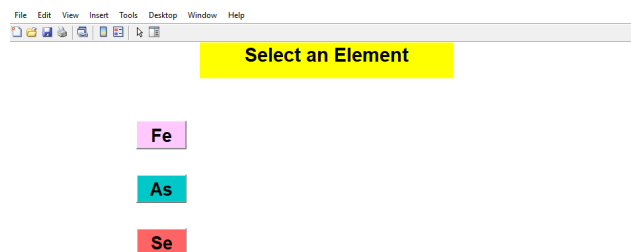

**Fig. 39**

Click on “open data files” and select the files for the standards in the folder Fe\_standards\_VP, then press “Plot”; do this for all Fe standards, except the “mixed valence” ones and you’ll obtain the scatter plot below (**Fig. 40**).

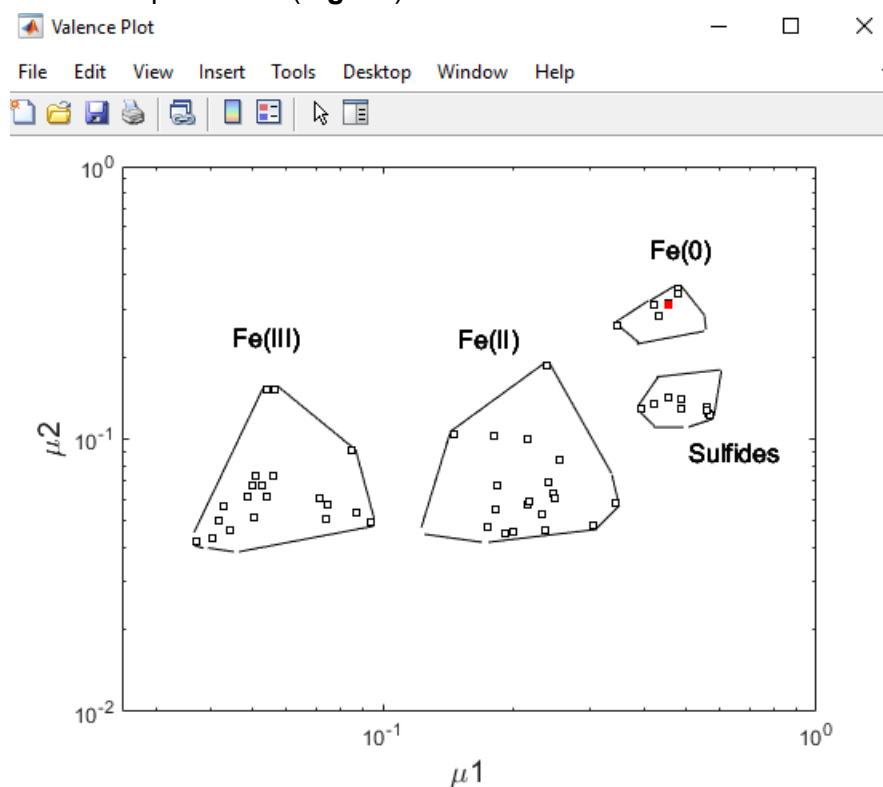

**Fig. 40**

The same holds true for the experimental spectrum that has been processed and saved as a txt file as explained in 3.2.2. Open the data file, select the desired plot options, then push “Plot” to see where it lies on the valence plot. You can choose several plot options listed under plot parameters (**Fig. 38**) circle, solid, red, font size 5 recommended. In the case it lands within the Fe (II) group, then it could be either a pure Fe(II) or a combination of Fe(III) with sulfides or metal. You’ll have to check your LSQ fitting results to assess as both go hand in hand. Note that the polygons are only indicative.

